

Address: JAMBU PERSHAD & SONS

6275/22 Nicholson Road, Ambala Cantt, Haryana, INDIA

Pin: 133001


Email: sales@japson.com japsonambala@yahoo.com

Website: www.japson.com **Phone:**

+91-171-4006897

Research Polarimeter

Product Image

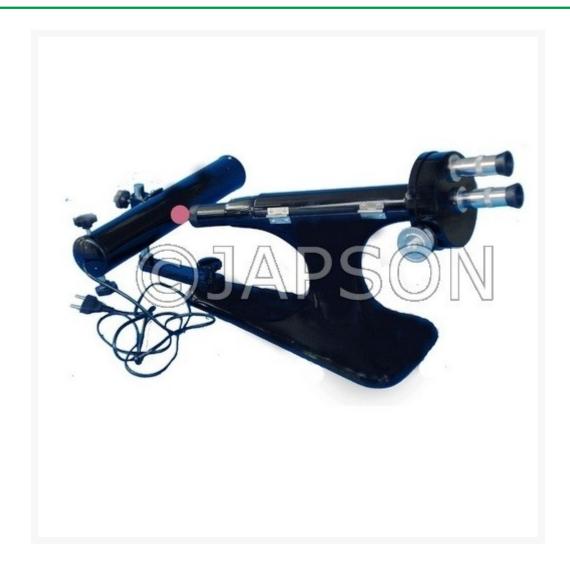

Description

Used in medical, pharmaceutical and chemical research laboratories. Consists of a imported half shade nicol polariser, analyser, 100 mm. diameter circular scale divided in degrees 0-360° x 1° and a sugar scale divided 30° to \pm 130° x 1° ISS, a divided micrometer drum is

given to read 0.05° (3 minutes) and 0.1° ISS. The scale is properly illuminated for taking readings. The body is very sturdy and separate sodium lamp stand is provided. Complete with polarimeter tube, bubble traps at the centre, 200 mm. Also available in 400mm model. Sodium lamp with stand and transformer.

Disclaimer

The Products details given on this page are indicative in nature and JAPSON reserves the right to change them without prior notice. Buyer is also requested to re-check the specifications and other features of product at the time of order as product development is a continuous process and minor modifications may be made to design based on latest availability, process and design.